

ZA180A Zero air

Ar zero

Código do produto: 30-1660 (230v) - UK | 30-1670 (230v) - EU | 30-1680 (110v) - US




Descrição

O gerador ZA180A fornece ar zero para aplicação como gás oxidante em detectores para GC, removendo hidrocarbonetos de uma fonte de ar comprimido e garantindo melhores resultados analíticos.



Generator Details

- **Tipos de gases:** Ar
- **Compressor integrado:** sim

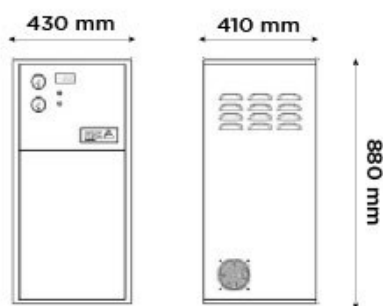

Características principais

- Fornece ar zero com nível de hidrocarbonetos inferior a 0,1 ppm
- Se necessário, pode funcionar 24 horas por dia, 7 dias por semana com excelente desempenho
- Um único gerador pode ser conectado a vários instrumentos
- Poucos itens consumíveis substituídos garantem baixo custo com manutenção ao longo de sua vida operacional
- Necessidade mínima de configuração
- Garantia de 12 meses

Especificações técnicas

	ZA180A
Tipo de gás	Ar
Fluxo máx. de gás:	18L/min
Teor de hidrocarboneto:	0.1ppm
Máx. pureza:	0.1ppm THC
Tempo de inicialização:	40 Minutes
Consumo de energia:	770 Watts
Voltagem:	110 / 230 ± 10% VAC
Frequência:	60 Hz
Corrente:	2/4 Amps
Temperatura máx. de operação:	30°C / 86°F
Certificações:	CE
Tamanho (AxLxP) mm	880 x 430 x 410 mm
Tamanho (AxLxP) pol.	34,6 x 16,9 x 16,1 inches
Peso do gerador	64kg / 140,8lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

